EN Electric Sonic toothbrush

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Rating: 100-240 V AC 50-60 Hz 0,9-1,3 W

Power supply: Li-ion accumulator 3.7 V

IMPORTANT SAFEGUARDS Read all instructions before use - Save these instructions

Basic safety precautions should be taken when using any electrical product

Toothbrushes can be used by children aged over 12 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.


To reduce the risk of electrocution: - Place and store the charging base at a safe distance from water and other liquids. - If the charging station falls into water, unplug immediately. Do not use it, but bring to authorized service center for fixing. - Do not reach for a product that has fallen into water. - Never use a charging station with a damaged cord or plug.

To reduce the risk of physical injury and to prevent product damage: - Do not use attachments other than those recommended by the manufacturer. - Do not plug the charging base into an outlet with a voltage other than the voltage specified at the bottom of the charging station. - Discontinue use if the product appears damaged in any way. - The charging station and toothbrush handle contain no user-serviceable parts. - Children should be supervised to ensure they do not play with the appliance

Medical Warnings - Consult your dental professional before use of this product, especially if you have had oral or gum surgery in the past 2 months.

- Discontinue use of this product and contact a physician/dentist if you experience discomfort or pain. - If you have a pacemaker or other implanted device, contact a physician or device manufacturer prior to use of this product.

DESCRIPTION OF PARTS

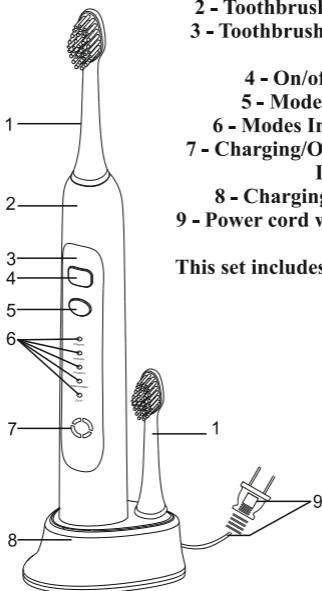1 - Toothbrush Heads
2 - Toothbrush Handle
3 - Toothbrush Control Panel

4 - On/off Button

5 - Mode Selector

6 - Modes Indicators

7 - Charging/Operating Indicator

8 - Charging Station

9 - Power cord with plug

This set includes 2 brush heads

CHARGING

Charge your toothbrush for 24 hours before first use.

This toothbrush will operate for up to 6 weeks on a single full charge.

- Place the toothbrush handle on the charging station. Then plug the power cord into an electrical outlet.

- Before each charging session, the lights will circle counterclockwise for 5 seconds to indicate the remaining battery level.


- The circle indicators consist of 4 lights, each equivalent to 25% of the total battery charge.

- 4 solid lights equal 100% of the remaining battery life.

4 rotating lights <25%

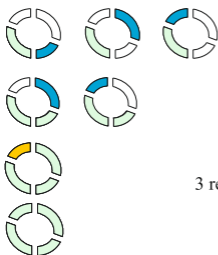3 rotating lights - lower left light is solid <50%

2 rotating lights - both lower lights are solid <75%

Upper left light is blinking
3 remaining lights are solid <95%

4 lights are solid 100%

TOOTHBRUSH USING

- Attach the toothbrush head to the toothbrush handle (Fig.1), leaving a small gap (about 0,5 mm.) between the two.
- Use your favorite brand of toothpaste.
- Guide the brush head slowly from one tooth to another in a small back-and-forth motion.
- Brush the gums as well.
- Do not scrub.

- A **2-minute timer with a 30-second pause** will remind you of when it is time to brush another section of your mouth.

This toothbrush can also be used by children over 12 years of age if they are supervised or given instructions on the safe use of the appliance. They must understand the hazards involved.

Cleaning and user maintenance shall not be performed by children without adult supervision.

BRUSHING MODES

CLEAN - up to 31,000 motions per minute, 3mm amplitude*.

Ultimate plaque removal.

WHITE - up to 41,000 motions per minute, 6mm amplitude*.

Extra power to help brighten and polish your teeth.

POLISH - up to 48,000 motions per minute, 3mm/4mm amplitude*.

MASSAGE - up to 31,000 motions per minute, 1mm to 3mm amplitude*.

Gentle mode for sensitive teeth and gums. Great for gum stimulation.

Keep pressing 'mode' button to switch between brushing modes.

DIY Mode offers five different vibration amplitudes for a personalized cleaning experience. To select the desired custom brushing mode, power on your toothbrush, then select 'custom' brushing mode by pressing the 'mode' button. Press and hold the 'on/off' button for 3 seconds until the 'custom' light blinks 3 times. Continue holding the 'on/off' button, then release at the preferred custom setting.

* **Of vibration**

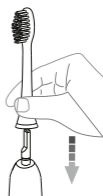
Fig.1

Toothbrush will power-on in the brushing mode used at the last brushing cycle.

AUTOTIMER This toothbrush is designed to reach the dentist-recommended brushing time of 2 minutes, with a pause auto-alert reminder every 30 seconds to change brushing areas.

Rotating light indicators will also inform you of remaining brushing times. After every 30 seconds of brushing, the corresponding individual light indicator will activate. After 2 minutes of brushing, all 4 lights will blink and the device will stop automatically.

30
seconds

60
seconds

90
seconds

2
minutes

TOOTHBRUSH CLEANING

- Clean your toothbrush regularly.
- Rinse the handle and the brush head thoroughly under running water.
- Remove the brush head from the handle and clean it separately (Fig.2). Wipe both dry before assembly and placing on base.

The shaft of the appliance becomes dirty from the toothpaste and water calc.

Clean the shaft with the cotton pad,

soaked with 40 % spirit solution once a week.(Fig.3)

Fig.2

Fig.3

CHARGING STATION CLEANING

- DO NOT clean the charging station by immersing it in running or still water.
- Before cleaning your charging station, make sure to unplug it from the electrical outlet.
- Use ONLY a dry or damp-dry cloth to wipe the surface of your charging station.
- Place the charging station onto a dry and safe surface until is completely dry.

NOTE! Manufacturer reserves the right to change technical characteristics of the device, with the aim of improving its quality, also change the parts and color of the appliance.

 **Disposal** Dispose the appliance in an environmentally friendly manner. It should not be disposed with the normal household waste. Please dispose it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

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